# NON SURGICAL INSTANT FACE-LIFT EMS DEVICE

## ENGLISH USER-MANUAL

Relieve fatigue  
Lifting & tightening the skin  
Diastolic & shrink pores  
Wrinkle & skin whitening

Convenient and Portable:  
Compact design for  
skincare on the go, It's like  
having a personal spa at  
your fingertips.

Apply to Face, eye, neck. Easy to carry, suitable for home and travel use!  
Invest in the future of your skin. Join countless others who have experienced  
the power of our EMS Face Massager. Get ready to embrace the natural beauty  
that's been waiting to shine through! Get your instant face-Lift effect now!

## INTRODUCTION

### 1. Acoustic vibrations beauty principle

Acoustic vibrations beauty principle... (The text in this block is partially obscured by a large watermark and is difficult to read fully.)

### 2. Hot and cool beauty principle

A Japanese professor at the university shows that repeatedly reaching the skin at a temperature of 42 °C helps to produce HSP protein, repair skin cells and strengthen skin elasticity to resist the formation of wrinkles. The warm 42°C treatment helps to dilute the pores for deep cleansing, while the cold treatment with a cold effect of up to 6°C instantly tightens pores and locks in nutrition. Alternating warm and frozen skin care, it stimulates the liquid circulation at the bottom of the skin and speeds up the metabolism of the cells. Besides effectively reducing the elimination of eye circles and simple eye lines, it also brightens skin. Even without makeup, the skin is also revealed. Combined with micro-vibration, it is beneficial to enhance the skin's firmness and relieve fatigue.

### 3. EMS beauty principle (Electrical Muscle Stimulation )

EMS Technology: Lift and tone your facial muscles for a youthful appearance. EMS is used to activate the skin cell, stimulating rejuvenation of collagen, to improve skin elasticity, soften wrinkles and leaves your skin looking smooth and healthy. It's like a workout for your face, helping you achieve that coveted, sculpted look without the gym and paying for expensive beauty salon treatments!

### Ice / COLD Effect:

1. Help skin firmness, maintain flexibility, prevent allergies and other problems;
2. Shrink pores, lock moisture and reduce wrinkles.
3. Can soothe muscles, relieve skin fatigue, soothe the skin sense of tension, reduce skin redness, pain, resistance to repeated growth of acne.

### Warm / HOT Effect:

By repeatedly exposing your skin to a temperature of 42°C, you can trigger the production of HSP. This extraordinary protein is renowned for its role in cellular health, promoting recovery and rejuvenation. The implications of this discovery are enormous and extend far beyond the realm of traditional skincare.


1. Warming, makes people feel warm and helps to relax muscles and relieve pain;
2. Stimulate the blood circulation acceleration, so the pores open, can help the skin to quickly import all kinds of skin care products absorb double as fast, make the skin more elastic and moist.
3. Putting a sponge to deep cleansing, better able to export the face toxins, blackheads, dirt.

### LED photon beauty principle

1. Red light: Wavelength 620nm(±10nm)18 LED's deeply reaches into the 630nm skin to promote blood circulation and the collagen regeneration enhance the skin elasticity, to make the skin to be smooth and moisture. Be suitable for the slack & wrinkly skin.
2. Blue light: Wavelength 470(±10nm) 18LEDs deeply reaches into the 480nm skin, to kill the bacteria, promote the blood circulation, can penetrate into 0.25mm under the cortex, destroy the bacteria living environment, anti-inflammatory, bactericidal, remove acne, acne scars and so on, effectively.

## NAME OF PARTS

- |                                 |                            |
|---------------------------------|----------------------------|
| ①. Hot and cool panel           | ②. LCD display             |
| ③. Charging jack                | ④. USB charging cable      |
| ⑤. Power & Mode switch          | ⑥. LED photon button       |
| ⑦. Level of EMS and Temperature | ⑧. Circle of makeup cotton |
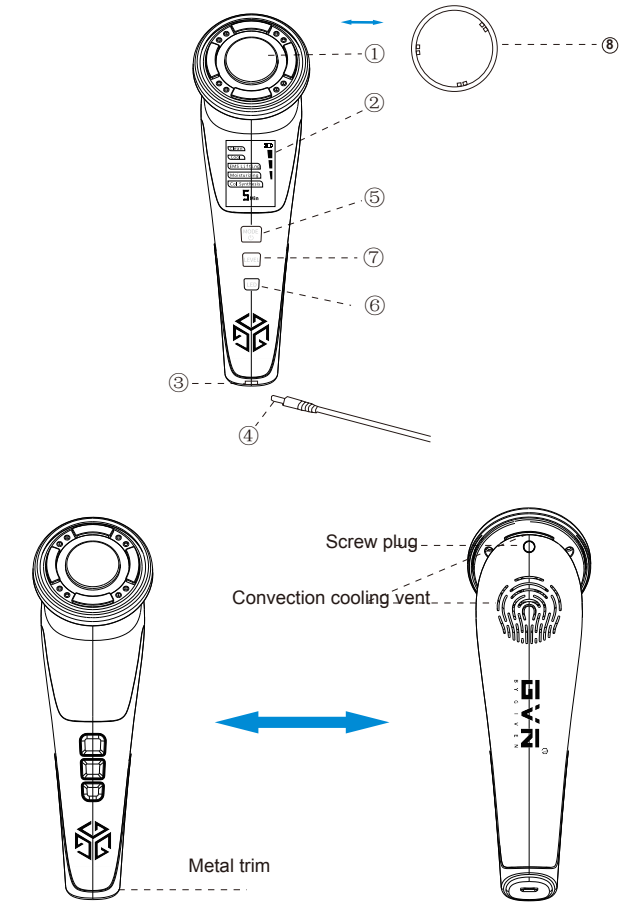
## SPECIFICATION

Name: Non Surgical Instant Face-Lift EMS device  
Input: DC 5V  
Vibration: 8000(±10)rpm  
Charging current: 500mA  
Battery capacity: 1000mAh  
Net weight: 145g  
Machine size: 175\*55\*67mm

Model: GVN 2021 V1  
Temperature: 6-42°C(±5°C)  
Charging time: 2.5Hours  
Power: MAX 8.5W  
Stand-by time: 45 days  
Use time: 5min/time or 2 min/time  
Box size: 215\*150\*80mm

## OPERATION INSTRUCTIONS IN DETAIL

1. Please fully charge the device before using it for the first time. The charging time is approximately 2.5 hours. Press and hold the power button for 5 seconds to start the device, with the default cleaning mode at level 1. The beauty device will then operate with vibrations and heat. Depending on your skin's needs, you can adjust the temperature by pressing button 7. You can adjust the LED Photon light by pressing button 6 in each mode (5 minutes of smart timing for mode 1).
2. Press the power button for the "Cooling Mode," and the beauty device will operate with vibrations and cooling. Depending on personal preferences, you can press button 7 to adjust the temperature (2 minutes of smart timing, after turning off, the fan will continue to work for an additional 2 minutes to cool down).
3. Press button 5 to switch to the "Hydrating Mode." The beauty device will vibrate and generate heat, along with EMS microcurrent to the skin. Depending on individual preferences, you can adjust the EMS intensity by pressing button 7 (5 minutes of smart timing for mode 4).
4. Press the power button 5 to switch to the "Collagen Synthesis Mode." The beauty device will vibrate and cool down while emitting EMS microcurrent to the skin.

According to individual requirements, you can adjust the EMS intensity by press key 7 (2min smart timing for mode 5).

### 6. Using Method as following pic and text:

	<b>Clean Mode:</b> Open pores and export pore dirt		<b>Collagen synthesis mode:</b> Lock skin nutrition, collagen synthesis
	<b>Cool Mode:</b> Shrink pores, tighten skin		<b>EMS Lifting Mode:</b> Stimulate muscle cells and repair wrinkles
	<b>Moisturizing Mode:</b> Promote nutrient absorption, lifting and tightening		

The "Cooling Mode" can effectively close the pores and maintain the skin's moisture balance and has a gentle cooling effect down to 6°C. Individuals who prioritize beauty are highly interested in the warm and cold treatment between 6-42°C. The "Heating Mode" can reach temperatures of up to 42°C.

EMS micro current stimulates the muscles, makes the skin more energetic and elastic, remove Crow's feet and wrinkles on the skin surface, and restores the skin's youthfulness, smoothness and tenderness.

After cleansing, apply skin care product such as GIVEN 24K B-TOX LAB Skin care products and GIVEN BTX-R Serum and activate the hydrating mode to promote the double absorption of skincare products! (In each mode, you can select red light/blue light depending on your skin's needs.)

After the skin care products were absorbed, turn on the collagen synthesis mode, cool to shrink pores, lock the nutrition. The EMS micro-current stimulates the skin cells to synthesize collagen and rejuvenate the muscles anti-wrinkle. (On each mode you can Select red/blue light according your skin needs.)

## SAFETY INSTRUCTIONS

1. Read all safety instruction before using.
2. Follow the operating instructions to use this device.
3. Do not pull down or remake the device, please consult the distributor for repairing.
4. Only use the device for the intended purpose as described in this manual.
5. Do not store or use the device in high-temperature and high-humidity places, such as in direct sunlight, near fire, in the bathroom, etc., to prevent fire or electric shock. Also, do not place it near heat-producing appliances.
6. Do not use or store the device in damp conditions to prevent short-circuits, electric shocks, and/or fires.
7. Keep it out of reach of children and physically disabled individuals, and only allow users under 18 years of age to use the device under close supervision.
8. When using the cooling function, the power consumption is relatively high and the fan has a delay for heat dissipation. It should be allowed to cool for about 10 minutes before it can be used again after using the cooling function.

### WARNING!

1. Please stop using when product abnormal heating, to avoid burns.
2. Please do not drop or clash the device.
3. Please apply serum constantly during use.
4. Please do not use it if the power cord is damaged.

## 5. Below conditions are not allowed to use! Heart disease, fever, communicable disease, cancer, open cuts, acute disease, hypotension, hypertension, install pacemaker, artificial heart or lung, during medical

## WARRANTY CARD

1. Free maintain for the breakdown caused by product quality problems.
2. One year guarantee from the date of purchase.
3. Free to change broken parts, freight is collected by customers.
4. If breakdown caused by the force majeure, free maintain won't be given and freight shall be collected by customers.
5. Please fill the Guarantee. If lost, won't be re-provided. guarantee clearly, without any revision.

Product Name	Item No.		
Customer Name	Phone		
Purchase Date	Means Of Purchase		
Address			
Records			
Date	Content	Results	Service Provider

Implementation of standards: GB 4706.1-2005

# NON SURGICAL INSTANT FACE-LIFT EMS DEVICE

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